



Handicare 2000

Twin rail curved stairlifts



The benefits of a stairlift from Handicare

If you have difficulty getting up and down stairs you are not alone. Every year thousands of people look for a solution to the problems posed by the difficulty they have using the staircase in their home.

Handicare's stairlifts have helped thousands of people overcome the problems posed by going up and down their stairs at home.


The decision to install a stairlift is sometimes a difficult one to take but when compared with the alternative of moving house, or living on one level of their house, many people have found it to be their preferred solution.

If you are thinking about installing a stairlift and your stairs turns corners, the Handicare 2000 offers you a safe and reliable solution. The stylish track system now gives you it's smoothest ever ride.

The Handicare 2000 stairlift suits most staircases and offers you all the basic functions that you would expect plus the option to upgrade to a powered folding footrest and/or powered seat swivel.

A tailor-made solution


Made to measure and perfectly finished it can be quickly installed at your home. After an informal visit we explain all the features of our stairlifts to you in detail.

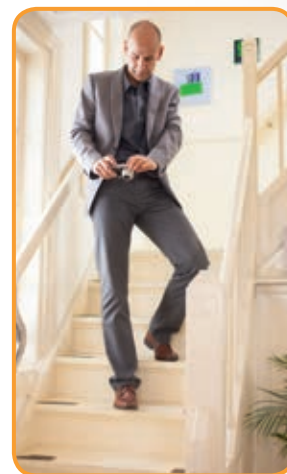

Style Seat

Each Handicare stairlift is designed to suit your individual staircase by using advanced PhotoSurvey3D technology, to make quick and accurate measurements. We have all the options available to ensure the best possible fit for your staircase.

Your Handicare 2000 stairlift can be fitted with a choice of seats and our Handicare supplier will be able to advise you which solution is most suited to your needs.

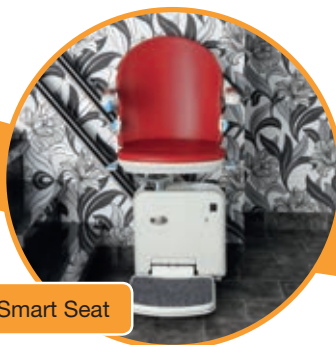

Smart Seat



Simplicity Seat



Standard Features



Remote control:

Your Handicare 2000 stairlift comes with two remote controls as standard, so you can have one upstairs and one downstairs. This means you can send the lift up and down the track without having to sit on it. This can also help free up space in your hallway or landing.



Safety as standard:

Your Handicare 2000 stairlift will come with a range of safety sensors that are designed to stop your stairlift if it meets any obstructions during travel. All Handicare stairlifts have a battery backup, which means it will continue to work in the event of a power cut.



Foldaway:

The stairlift folds away, leaving the staircase free for other users. The tailor made track is designed to maximise space on the stairway.





Easy to use:

The operating toggle is very easy to use, simply hold it towards the direction you wish to travel and your stairlift will move smoothly on its way.



Seat belt:

For added peace of mind, your Handicare 2000 is fitted with an easy to use retractable seat belt. If you prefer, you can choose the Velcro seat belt solution.



Diagnostic display:

Your Handicare 2000 stairlift has a clear diagnostic display, so that you can monitor the status of your stairlift, for example if a safety sensor has been activated.



Comfort features

Automatic hinge track

Some homes will have a doorway close to the bottom of the stairs. An automatically hinged track is available for these situations and means that the track will be folded away from the doorway to give clear access. When operating the stairlift from the toggle, the hinge will lift and lower automatically.



Powered seat swivel

When you arrive at the top of the stairs, you will need to swivel your seat round, so that you exit onto the landing safely. If you choose the powered option, the swivel is automatically operated by continuing to hold the toggle when you arrive at the top of the staircase.




Weight capacity

The standard maximum weight limit on the Handicare 2000 is 250 lbs. This weight limit can be increased to 300 lbs. by having the heavy duty kit fitted on your stairlift.

Powered footrest

The footrest and seat can be folded up to leave the staircase as clear as possible for other users to walk past. If you choose the powered footrest option it will remove the need for you to bend down to fold or unfold the footrest.

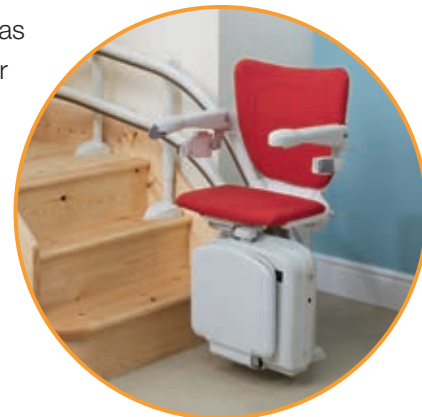

Stairlift colors

Fitted to the stairs and not the wall, it is unlikely you will need to redecorate following the installation of your Handicare 2000 stairlift. With a choice of upholstery and track colors, your stairlift will compliment most color schemes. For an additional fee, the rail can be finished in any color to suit your décor.

Perch seat

If you have a narrow staircase or if restricted movement in the knee or hip joint makes sitting painful, a perch seat may be the solution.

A perch seat means that the user is almost in a standing position when using the lift. There is a shortened seat and retractable seat belt which provide a stable support area for users to 'perch' on. The upholstery is available in a choice of six colors.




Handicare 2000 series: at a glance

Opzioni sedile	2000 Smart	2000 Simplicity	2000 Style	
Choice of Upholstery	✓	-	✓	A selection of six or three colors means you can choose upholstery to match your décor.
Powered Swivel	○	-	○	This makes turning the stairlift seat to get in and out on the landing as easy as possible.
Retractable Seat Belt	✓	✓	✓	This retractable belt can be clipped together on one side or centrally.
Velcro Seat Belt	○	○	○	This black belt has velcro on either end to hold it firmly in place.
Perch Seat	○	-	-	Ideal for users who have difficulty bending at the hip or knee joint.
Footrest Options	2000 Smart	2000 Simplicity	2000 Style	
Manual	✓	✓	✓	The footrest is folded and unfolded manually.
Powered from Seat	○	○	○	Folds the footrest away when the seat is folded up.
Powered from Arm	○	-	-	The footrest can be lifted and lowered by a switch under the arm whilst the user is seated.
Track Features	2000 Smart	2000 Simplicity	2000 Style	
Travels around corners and bends	✓	✓	✓	For a wide range of staircases, from ones that turn a corner to multi-flight or spiral.
Multiple charge points	○	○	○	Use your remote control to park on charge points.
Track available in any color	○	○	○	Customise your curved stairlift track to match your décor.
Zero intrusion into landing	○	○	○	For use when a doorway close to the top of the stairs would be obstructed by the track.
Powered hinge	○	○	○	Automatically moves the track while the stairlift moves leaving doorways unobstructed.
Weight capacity	2000 Smart	2000 Simplicity	2000 Style	
Standard Maximum Weight	250 lbs	250 lbs	250 lbs	
Max weight with Heavy Duty Kit	300 lbs	300 lbs	300 lbs	
Standards compliance	2000 Smart	2000 Simplicity	2000 Style	
EN 81-40: Stairlifts	✓	✓	✓	


✓ Available as standard ○ Optional - Not available

Smart Seat: Upholstery Color Choices

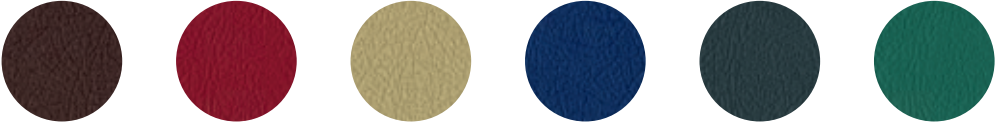

COCOA RUBY SAND SAPPHIRE SLATE JADE

Simplicity Seat: Upholstery Color




SAND

Style Seat: Upholstery Color Choices



CHERRY RED VANILLA BEIGE IVORY WHITE

Track Color Choices:



LIGHT GREY BROWN IVORY WHITE

Handicare also offers special orders so the track can be finished in any color to suit your décor.


Please be aware that actual colors may vary slightly from those shown.



Making everyday life easier

Handicare Vision: a 3D preview of a stairlift on your stairs

When it comes to stairlifts we want you to make the best choice possible. You can have a 3D preview of Handicare's stairlift range in the comfort of your own home, and see how each one would look in your hallway too, thanks to Handicare's Vision app. Vision is a special application we have developed to show you what a stairlift would look like on your stairs, and to help you make the right decision.



Quality and safety

All our new products are delivered with a comprehensive warranty.

A Handicare stairlift stands for high quality: quality in terms of safety, innovation, comfort and design. Furthermore, all our products are tested and approved by worldwide independent inspection services. This gives you the assurance that all our products meet the highest standards.



Experience

Convenience serves people... and lifting equipment of Handicare has been doing just that for more than 125 years.

With technical perfection, intelligent solutions and safe applications, Handicare's innovations have provided not only greater comfort, but also more attractive stairlifts... in your style, your taste.

Top quality in all respects make the Handicare stairlift a valuable asset for independence around the home. Many satisfied users experience that every day.



Handicare
2201 Hangar Place, Suite 200
Allentown, PA 18109
<http://handicarestairlifts-us.com>
customerservice.usa@handicare.com

Local 610-266-5260
Toll free 866-276-5438
Fax 610-266-5266



Handicare Freecurve

Single Rail Curved System



handicare

The benefits of a stairlift from Handicare

Handicare has produced and installed lifts for more than 125 years and has many satisfied and happy users worldwide. Our innovations have created some of the most comfortable, safe and quiet stairlifts. A curved stairlift from Handicare is not only convenient, but also looks great and fits almost any home with its variety of available seats and rail colors. Top quality in all respects, a Handicare stairlift will be a convenient and valuable asset.


You will have independence with contemporary design and maximum safety for a price that will certainly appeal to you. With a one-off investment you will regain your mobility and be able to stay in your own home.

A Handicare stairlift will provide you comfort, safety and confidence to again enjoy the place you love.


A perfect solution for every staircase

Handicare has a solution for every staircase. All curved stairlifts share the same typical features: easy to operate, whisper silent and functional design with the single tubular rail.

With our different types of stairlifts we are able to place a Handicare stairlift in almost every type of house and provide a safe and comfortable ride up even the most complicated of staircases. The compact design allows plenty of space for other stair users.

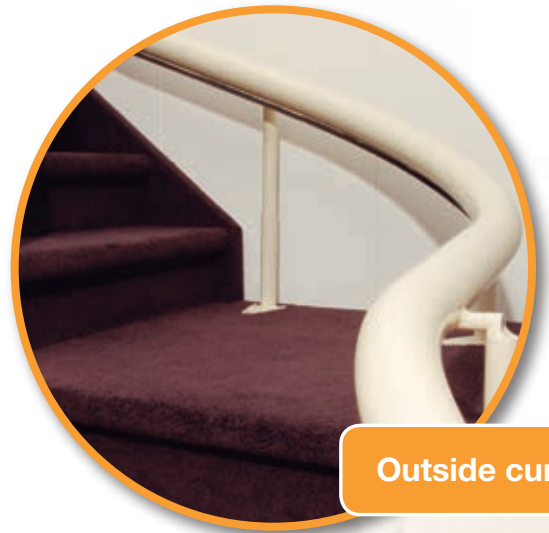


Outside curve

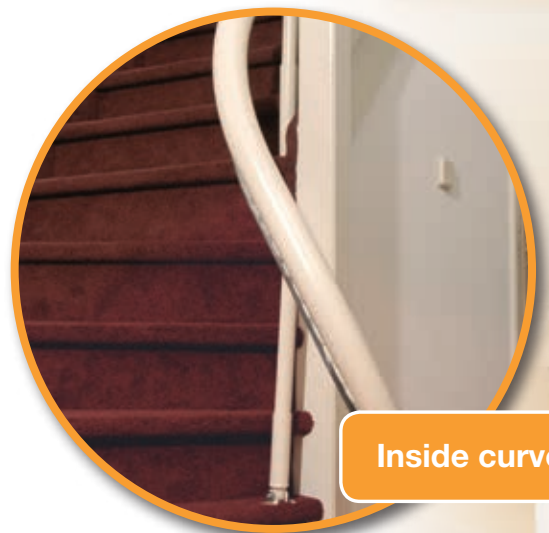

Inside curve



Multiple floors

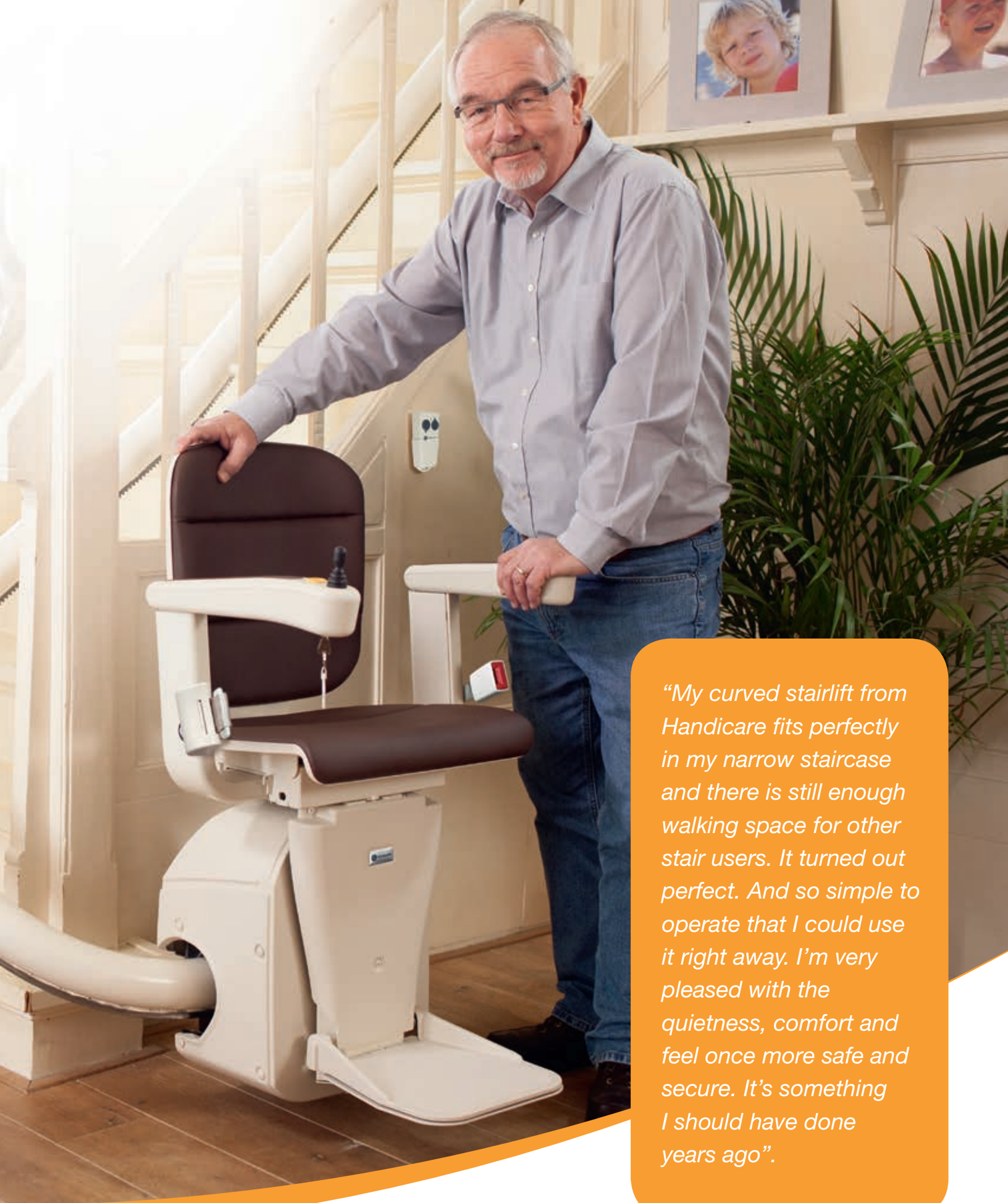

“My curved stairlift from Handicare fits perfectly in my narrow staircase and there is still enough walking space for other stair users. It turned out perfect. And so simple to operate that I could use it right away. I’m very pleased with the quietness, comfort and feel once more safe and secure. It’s something I should have done years ago”.

A seat for every interior

Handicare offers you the choice of different stairlift seats, each of them available in a range of colors. There's plenty of choice, in terms of both comfort and style to suit that of your home. Because convenience can look good too.



Classic

The best-selling Classic seat has an attractive leather grain effect finish. Made from durable, non-slip polyurethane, it is hard wearing and easy to keep clean.



Alliance

The cream upholstery finish accompanied with the dark wood seat frame gives the new Alliance seat it's luxurious and stylish look. Truly a seat that matches every interior.

Upholstery Color Choices*



ANTHRACITE BORDEAUX
RED BEIGE

Manual swivel

Push down lever allows seat to be swiveled enabling safe mounting and dismounting.



Powered swivel**


The seat will move automatically, operated by the joystick.

Seat features




Remote control

With light pressure on the call and send buttons the lift is always within arm's reach.

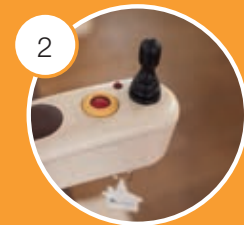

Joystick

This easy to use control can be positioned on the left or right arm.

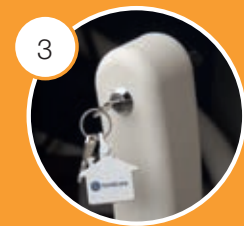

Key switch

The key can be removed to prevent unauthorized use.



Safety belt

Fitted as standard for your protection and security.



Footrest linkage**

Simultaneously folds up with seat.



Elegance

The Elegance seat is lightly padded and offers a harmonious combination of stylish design and sophisticated comforts.

Upholstery Color Choices*



BORDEAUX RED



COCOA BROWN

Space saving design

The arm rests and foot support can be folded away easily, maximizing the clear space on the stairs.



Safety sensors

Your Handicare stairlift will come with a range of safety sensors that are designed to stop your stairlift if it meets any obstruction as you travel up or down the stairs.

* Please be aware that actual colors may vary slightly from those shown
** Optional



A tailor-made solution


Perfectly finished, 100% made to measure and yet rapidly installed in your home. Your dealer, after an informal visit, will explain all the features of our stairlifts to you in detail.



Each Handicare stairlift is designed to suit your individual staircase. We use advanced PhotoSurvey technology to make accurate measurements and we have all the options available to ensure the best possible fit for your staircase.

Handicare Vision: a 3D preview of a stairlift on your stairs

When it comes to stairlifts we want you to make the best choice possible. Therefore, we are delighted to offer you additional help to make choosing the right stairlift a lot easier. Vision is a special tablet application developed by Handicare. With this application your stairlift advisor can show you how your Handicare stairlift will look on your staircase. Vision uses augmented reality to show you 3D images of various seats, colors and options from our assortment. With the 3D preview you can see and choose your Handicare stairlift in the comfort of your own home.

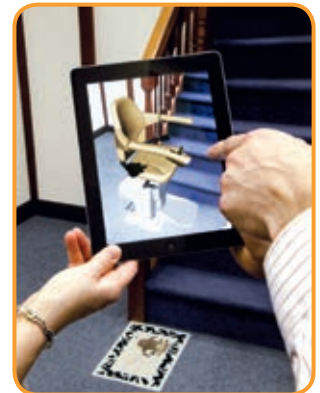


Solutions for narrow stairs

Where a stairwell is particularly narrow, or a user has limited knee flexion, you can fit a curved stairlift with a 45° swiveled seat or with the 70° swiveled seat of the Turn & Go system. This enables curved stairlifts to be installed on stairs of only 27 1/8" wide and could avoid the inconvenience of major building adaptations, or save on the cost of a 'through the floor' lift.




Automatic folding hinge

To avoid any obstruction at the bottom of the stairs, e.g. where there is a doorway, an automatic hinged rail with safety system is available.

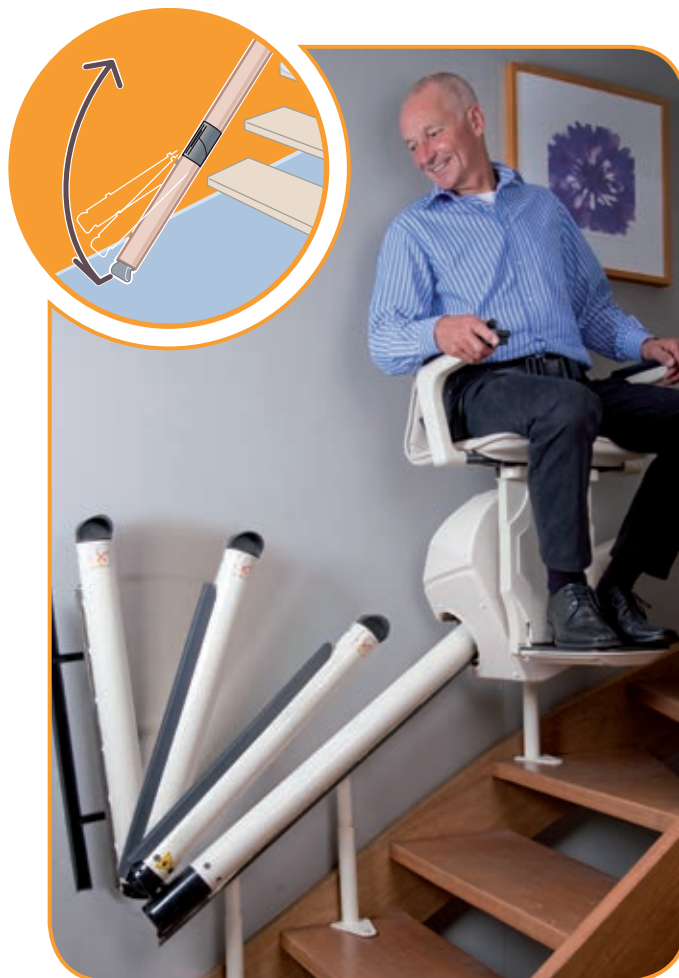

A rail color that matches

Handicare's curved rail matches every interior and is available in cream-white, brown and grey. For an additional fee the rail can be finished in any color to suit your décor.

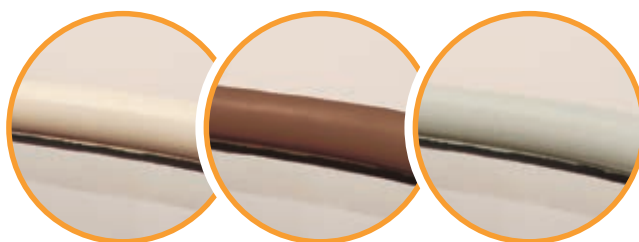

Making everyday life easier

Experience

Convenience serves people and the lifting equipment of Handicare has been doing just that for more than 125 years.

With technical perfection, intelligent solutions and safe applications, Handicare's innovations have provided not only greater comfort, but also more attractive stairlifts to compliment your home.

Top quality in all respects make the Handicare stairlift a valuable asset for independence around the home. Thousands of satisfied users experience this every day.

Quality and safety

All our new products are delivered with a 2 year warranty.

A Handicare stairlift stands for high quality: quality in terms of safety, innovation, comfort and design. Furthermore, all our products are tested and approved by worldwide independent inspection services. This gives you the assurance that all our products meet the highest standards in terms of quality and safety.



Handicare
2201 Hangar Place, Suite 200
Allentown, PA 18109
www.handicare.com/us
customerservice.usa@handicare.com

Local 610-266-5260
Toll free 866-276-5438
Fax 610-266-5266

Please be aware that actual colors may vary slightly from those shown. All information is provided in good faith but does not constitute a guarantee. Handicare reserves the right to alter designs and specifications without prior notice.