




"The ideal match for  
my active lifestyle"

Sarah Williams

## The perfect fit for your lifestyle

---

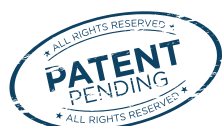

reddot award 2015  
winner




Yes, you can.®



## The perfect fit for your lifestyle

The *Invacare MyOn<sup>HC</sup>* has been designed for the more active user looking for an affordable chair with good drive performance and the ability to fold easily for handling, storage and transportation purposes. With its two patents pending and lightweight construction, the new *MyOn<sup>HC</sup>* gives the user a chair they can rely on. Key features include:

- **Unique folding mechanism**
- **Stability of a rigid chair**
- **Quick step-less adjustability**
- **Highly configurable**

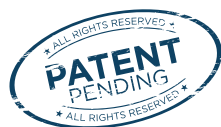

reddot award 2015  
winner



**Yes, you can.®**

## ARMRESTS



1



2



3

## NEW SWING AWAY AND REMOVABLE ARMREST



4



5



6



7



8

### ARMRESTS

#### 1 Mudguard

The stylish mudguards are lightweight keeping the overall weight of the wheelchair to a minimum.

#### 2 Conventional height adjustable armrest

Depth adjustable by sliding on the seat rail.

#### 3 Hemi-armrest

This optional armrest provides additional comfort for the user with its height, depth and angle adjustability.

### NEW SWING AWAY & REMOVABLE ARMREST

#### 4 New, swing away and removable armrest

Accessible by the user and swings out to assist with transfers.

#### 5 Height adjustment

Simple push button providing easy adjustment.

#### 6 Release lever

Even with limited hand function, the armrest can easily be released and flipped backwards.

#### 7 Flip back position

A simple screw setting allows the position of the armrest to be fixed for easy access when flipped backwards.

#### 8 Reversible armpad

The position of the armrest can be easily reversed providing a short and long armpad in one.



### Unique folding mechanism

The innovative design of the chair allows the frame to **greatly reduce in size** and the backrest to fold on top of the seat, making the chair **compact**

No additional handling is required to fold the chair, making it effortless for **storage and transportation**



### Stability of a rigid chair

The patent pending 'H-Shaped' horizontal folding mechanism provides an **improved level of stability** compared to other folding chairs on the market

The self-locking function ensures a secure seat position and makes the chair almost effortless to propel, **improving drive performance**



### Adjust the chair whilst the user remains seated

Our unique, patent pending technology allows backrest angle adjustments to be made whilst the **user remains in the chair**

The step-less adjustability makes **assessments and fittings both quick and easy** to ensure the perfect fit for the user

The built in scale allows for a **more accurate backrest angle** to be set



## Make it personal

Giving the user the freedom to design the chair around their personality, there are a variety of colours available to choose from

The **MyOn<sup>HC</sup>** can also be **configured depending on the needs of the user**, starting from selecting a more active look with the fixed frame, or the swing away version for even easier folding and transportation

The next two pages detail the range of options available on the **MyOn<sup>HC</sup>**

The **MyOn<sup>HC</sup>** is also **cross compatible** with other Invacare products including the Alber power add on range.



## Award winning

Already an **award winning** chair, the Invacare **MyOn<sup>HC</sup>** has won a Red Dot Award for Product Design

The Red Dot Award is not only the largest product competition in the world, but also one of the most prestigious design awards




**reddot** award 2015  
winner



## CONFIGURABILITY



1

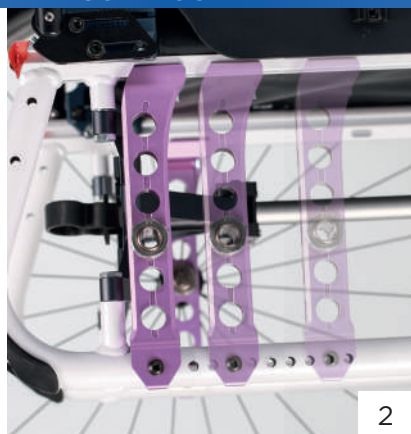

2



3

## LEGREST OPTIONS



4



5



6

## ACCESSORIES



7



8



9

### CONFIGURABILITY

#### 1 Tension and depth adjustable seat

Adjustments can be made via easy to use straps.

#### 2 Centre of gravity

100mm centre of gravity adjustments can be made in 25mm increments which can even be fine-tuned to 12mm.

#### 3 Camber

Standard adjustments of 0-2° available with the option of 3 & 4° adjustments if required.

### LEGREST OPTIONS

#### 4 Swing away legrests

Allows the leg rests to swing either inwards or outwards.

#### 5 Fixed front

An 80° fixed front frame is available for the more independent users.

#### 6 Selection of legrest options

Choose from 70, 80 or 90° leg rests, elevating legrests, or an amputee support.

### ACCESSORIES

#### 7 Backrest cover

Includes a backrest pocket and reflective stripe for better visibility.

#### 8 Active anti-tipper

Can be operated by both the user or the caregiver for additional support.

#### 9 Front wheels

A selection of wheels are available from 4" to 7", including stylish aluminium soft roll options.

## TECHNICAL DATA

Invacare <b>MyOn<sup>HC</sup></b>		12", 13", 14", 15", 16", 17", 18", 19", 20" and 21"		14", 15", 16", 17" 18", 19" and 20"		0 - 9 °		Front 14"-21" in steps of 1" Rear: 14" to 20" in steps of 1"		Fixed 90° Or stepless 90° -15°/+15°		Fixed 12", 13", 14", 15", 16", 17", 18", 19" or 20" Adjustable from 14" -16" or 16"-19.5"		Fixed front 80° Swing away legrest 70° / 80° & 90° Elevating legrest 0-80°
--------------------------------------	--	---	--	--	--	---------	--	---	--	---	--	---	--	--

Invacare <b>MyOn<sup>HC</sup></b>		Swing away frame: 8"-20" Fixed frame: 12"-21"		6" - 10.5" 8" - 11" 10" - 15"		SW +7" To SW+ 10"		12.5" - 15"		33" - 47.5"		26.5" - 34"		27.5" - 43.5"
--------------------------------------	--	--	--	-------------------------------------	--	----------------------	--	-------------	--	-------------	--	-------------	--	---------------

Invacare <b>MyOn<sup>HC</sup></b>		7°		41.5" - 54.5"		220 lbs & 290 lbs		From 27 lbs		From 18 lbs
-----------------------------------	--	----	--	---------------	--	-------------------	--	-------------	--	-------------



This product was tested and measured according to RESNA WC\_1:2009. Dimensions may vary by 15 mm/2°. For more information about this wheelchair, please visit [www.invacare.ca](http://www.invacare.ca). For information about typical user description, intended use, adjustments and transportation, please refer to the user manual.

## FRAME COLOURS

Deep Red	Pop Blue	Anthracite Black	Pearl White	Extreme Green	Hypnotic Purple
Pearl Grey	Matte Black	Intense Yellow	Vibrant Orange		

## UPHOLSTERY



Black Nylon /  
Grey Mesh

## COLOUR ACCENTS

Silver	Blue	Purple

Please notice colours may vary slightly from those displayed above.

**Invacare Canada**  
[www.invacare.ca](http://www.invacare.ca)

Canada  
570 Matheson Blvd. E.,  
Unit 8  
Mississauga, Ontario  
L4Z 4G4 Canada  
(800) 668-5324

©2018  
Invacare Corporation.  
All rights reserved.  
Trademarks are  
identified by the  
symbols <sup>™</sup> and <sup>®</sup>. All  
trademarks are owned  
by or licensed to  
Invacare Corporation  
unless otherwise  
noted. Specifications  
are subject to change  
without notification.

Form No. 18-375C



**Yes, you can.<sup>®</sup>**